



**EBERLE**

Hip-perfect burr

# Hip-perfect burr

by Prof. Dr. med. Stefan Landgraeber



## Single use

shaft diameter: 4.2 mm  
burr diameter: 5.0 mm  
effective length: 210 mm  
bend: 20°



## Eberle

article no. 9143642H



## Arthrex and Conmed/Linvatec


article no. 023642HS

In the arthroscopic treatment of femoroacetabular impingement of the Cam type, precise resection of the bony protrusions on the head-neck junction near the joint is a particular challenge.

In order to avoid under or over-correction, a good overview is needed. While this is normally ensured in the ventral portion of the femoral neck, the overview may be limited during milling at the lateral and especially in posterolateral regions.

The design of the Hip-perfect burr enables a good view even of these hard-to-reach places due to its special curvature and narrow design. The especially machined ball head enables reliable milling in nearly all positions.

**EBERLE**

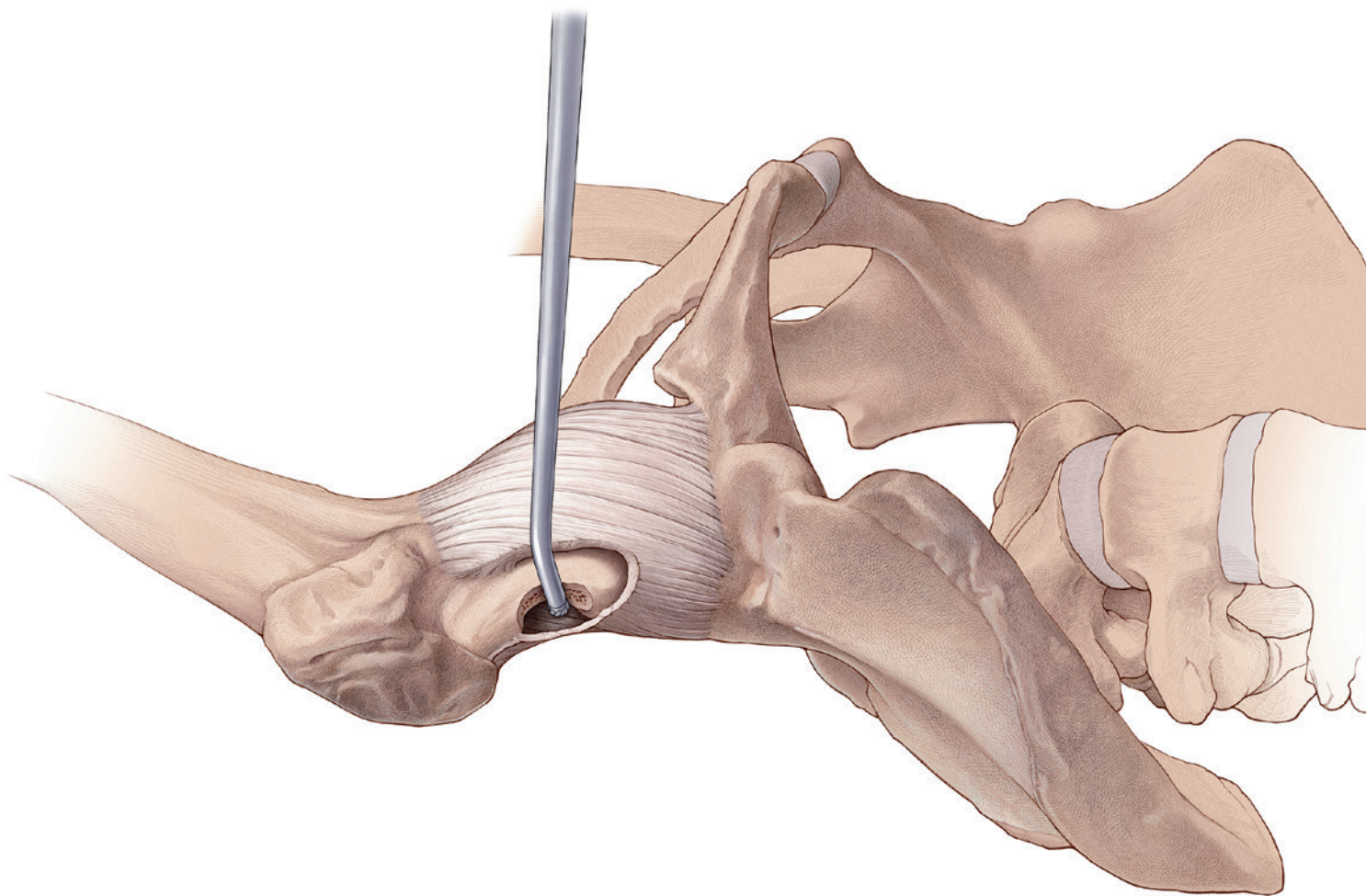

**EBERLE**

**EBERLE GmbH & Co. KG**

Dachsteinstraße 2  
75449 Wurmberg  
Germany

Phone +49 (0) 7044 9611-288

[sales@eberle-med.com](mailto:sales@eberle-med.com)  
[www.eberle-med.com](http://www.eberle-med.com)

---

